Global Diagnostics Supplier





UBM - F_x-19ET Auto Hematology Analyzer

Better user–experience with single chamber (21 parameters and 3 histograms)



Specifications

Parameters: WBC, LYM#, MID#, GRAN#, LYM%, MID%, GRAN%, RBC, HGB, HCT, MCV, MCH, MCHC,RDW-CV,

RDW-SD, PLT, MPV, PDW, PCT,P-LCR,P-LCC

3-part differentiation of WBC, 21 parameters and 3 histograms

Principles of operation:WBC/RBC/PLT: Electrical Impedance

HGB: Photoelectric colorimetry

Aspiration Volume: Whole blood 10ì L, Pre-diluted 20ì L

Aperture Diameter:WBC/RBC/PLT 80i m

Throughput:30 samples per hour

Data storage: More than 100,000 sample results with histograms

Alarms: Error messages

Display: 10.4 inches color LCD, resolution 640 ×480

Control and calibration:L-J,X,X-R,X-B four control modes,manual

and auto-Calibration

Power requirement: AC 100V --- 240V, 50/60Hz, 130VA-180VA

Input/Output: Support RS232, USB port,

Shortcut keys with standard mouse and keyboard


Internal 57mm thermal printer,

Parallel Port and USB for external printer





System check

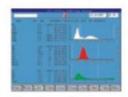

Counting results



System Maintenance

100		STREET, 17 hard	WCA5
NACON.			
	-		
	_	NOT WELL	
-		COLUMN TRACE	
-	-	Mr. hallan	
_	-	Section 1995	
-		to accompany	
		Miller and Control	
		Student region.	
-		Other Review	
-		Control of the Contro	
		100 700	
		100 001	
		The second	

System Log

Characteristics

- Large color LCD display & Linux operation system
- Support USB data backup & system upgrade
- Excellent data management
- High efficient self checking system, less maintenance
- Support LIS and HIS with HL7 protocol
- Collaborator of BIO-RADcontrol, target value available

Unitron Bio-Medicals

AN ISO 9001: 2015 CERTIFIED COMPANY

Mfs., Sales & Services Instruments & Diagnostics Reagents
"Shree Durga", # 464/1, 12th Main Road, 17th 'A' Cross,
J.P. Nagar 5th Phase, Bangalore-570 078. INDIA.

Phone: +91-80-25187600

E-mail: info@unitronbiomedicals.com : ubm06@rediffmail.com Web: www.unitronbiomedicals.com

Customer Care: 94490 09531